



Viral DNA and RNA Extraction Kit

User Guide

Version 8.0



In-Vitro Diagnostics / For use with Automatic nucleic acid extractor compatible with
Viral DNA and RNA Extraction Kit



T014H T016H T060H




Xi'an Tianlong Science and Technology Co., Ltd.


No.4266, Shanglin Road, Weiyang District, Xi'an, 710021, Shaanxi, P.R. China




SUNGO Europe B.V.

Fascinatio Boulevard 522, Unit 1.7, 2909VA Capelle aan den IJssel, The Netherlands

 T014H
64

 T060H
32

 T016H
20

Contents

Intended Use	1
Product Performance Indicators.....	1
Special Notes.....	1
Testing Principle.....	1
Content of the Kit.....	2
Materials Required but not Provided.....	2
Warnings and Precautions.....	2
Precautions for Safe Handling.....	4
Reagent Storage and Handling.....	4
Sample Handling and Storage.....	4
Operation Guide.....	4
1. Automated Extraction Process.....	4
2. Operation Steps of Automated Extraction.....	4
2.1 Automatic Nucleic Acid Extractor (model: Libex).....	4
2.2 Automatic Nucleic Acid Extractor (model: GeneRotex 96).....	6
Troubleshooting Guide.....	8
Quality Control.....	9
Limitations of Test Methods.....	9
Safety Symbols and Signs.....	9
Contact Information.....	10

Kit Version	8.0		
Changes	Address of Manufacturer Address of EU Representative Chapter "Intended Use" Chapter "Content of the Kit" Chapter "Materials Required but not Provided" Chapter "Warnings and Precautions" Chapter "Sample Handling and Storage" Chapter "Safety Symbols and Signs" Small lexical corrections.	Additions	Chapter "Precautions for Safe Handling"

Intended Use

The ***Viral DNA and RNA Extraction Kit*** is designed to rapidly extract viral DNA and RNA from swab samples. The extracted viral DNA and RNA are of high purity and stability and can be used in a variety of routine operations, including enzyme digestion, Polymerase Chain Reaction (PCR), DNA library constructions, Southern hybridization and blotting and other experiments.

The ***Viral DNA and RNA Extraction Kit*** is intended to be used by professionals, such as biotechnologists, microbiologists, clinical technicians, and physicians who are trained in molecular and biological techniques.

Product Performance Indicators

The ***Viral DNA and RNA Extraction kit*** can extract nucleic acids from swab samples in high-efficiency, particularly from low-copy complex samples. Both the intra and inter-batch variations of kit are less than 5%.

Special Notes

The ***Viral DNA and RNA Extraction Kit*** is worked with TIANLONG® automatic nucleic acid extractors (Libex and GeneRotex 96) that have been disinfected by UV light before use. After an experiment, wipe the inside of the extractor with 75% ethanol and disinfect it with UV light for 15 mins. An automatic nucleic acid extractor automates the entire purification process and can process 1-96 samples in a single run.

The kit is used to extract viral DNA and RNA targets. To avoid RNA degradation by RNase during operation, use exclusive-use utensils and sample injectors, and use disposable centrifuge tubes and tips processed by autoclave before using. The operator should wear powder-free gloves and a mask and a protective cover all.

The kit has magnetic beads with a unique separation function and a unique buffer system to extract, isolate and purify high-quality nucleic acids from a variety of liquid samples such as swab samples.

Magnetic beads enable the purification of high-quality nucleic acids that are free of protein, nuclease, and other impurities. Purified nucleic acids can be widely used in variety of routine operations, including downstream experiments such as enzyme digestion, Polymerase Chain Reaction (PCR), DNA library construction, Southern hybridization and blotting.

Please carefully read the manual of instructions before attempting to install or use the product for the first time. To consider all possible consequences of incorrect operation or non-recommended functions, pay special attention to the possible consequences.

Testing Principle

The ***Viral DNA and RNA Extraction Kit*** is worked with TIANLONG® automatic nucleic acid extractors (Libex, GeneRotex 96). During the nucleic acid extraction process, magnetic beads are adsorbed, transferred and released using special magnetic rods based on the principle of C adsorption. The extraction process enables the conduction of nucleic acid extraction and final adsorption of highly pure nucleic acids with the transfer of magnetic beads and nucleic acids.

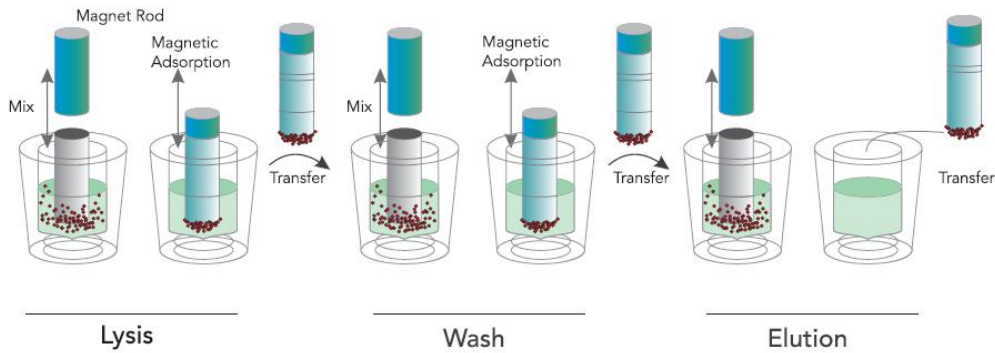


Figure 1. Schematic Diagram of Automatic Nucleic Acid Extractor

An automatic nucleic acid extractor performs the following steps on a sample containing magnetic particles:

A magnetic rod protected by the mixing sleeve inserts into a well which contains sample. The mixing sleeve stirs rapidly and repeatedly in the liquid to ensure complete mixing of the liquid and magnetic beads. After cell lysis, nucleic acid adsorption, washing and elution, highly pure nucleic acid is obtained.

GeneRotex 96 equipped with an array of 96 magnetic rods, allowing it to process up to 96 samples simultaneously.

Content of the Kit

Short Code		T014H	T016H	T060H
		Name of Component		
REAG1	Size	64 T/Box	20 T/Box (DT6)	32 T/Box
	Component	Pre-filled 96-deep well plate	Pre-filled 6 strip tube	Pre-filled 96-deep well plate
	Quantity	4	20	4
	Component Specification	16 Tests	1 Test	8 Tests
Instructions for Use		1 Copy	1 Copy	1 Copy

Materials Required but not Provided

When working in a laboratory, make sure to wear a proper lab coat, powder-free disposable gloves and protective goggles. For more information, please consult the Safety Data Sheet (SDS) available from the product supplier.

- Pipettor: 200 µL or 1000 µL
- Tip: 200 µL or 1000 µL
- Vortex mixer
- Sample holder
- 75% ethanol
- Single kit docking (matched with T016H (6 strip tube), could be purchased from Tianlong)
- Extractor

Warnings and Precautions

Please be sure to read the precautions before using the kit.

The kit is used to extract viral DNA and RNA targets. To avoid RNA degradation by RNase during operation, use exclusive-use utensils and sample injectors, and use disposable centrifuge tubes and tips processed

by autoclave before using. The operator (researcher or clinical expert) should wear powder-free gloves and a mask.

Please read the manual carefully before using the kit, and strictly follow the manual throughout the operation. The clinical samples should be collected on a clean bench or in a bio-safety cabin.

Before using TIANLONG® automatic nucleic acid extractors (Libex and GeneRotex 96), they must be disinfected by UV light. After an experiment, wipe the inside of the extractor with 75% ethanol and disinfect it with UV light for 15 mins.

Due to the possibility of residual magnetic beads in the eluate following extraction, every possible effort should be made to avoid suctioning of any magnetic beads during eluate absorption.

Do not mix reagents from different batches and use the kit within expiry date.


Dispose of all samples and reagent materials used in an experiment, thoroughly clean and disinfect the experimental work bench.

The **Viral DNA and RNA kit** is intended for in vitro diagnosis use.

When using kit, always wear a suitable lab coat, disposable gloves, and protective goggles. For more information, please consult the appropriate Material Safety Data Sheets (MSDSs). These documents are available online in a convenient and compact PDF form at <https://www.medtl.net/resources/download/catalogue-all/catalogue>, where the operator can find, view and print the appropriate MSDSs.

⚠ Caution: Do not add any bleach or acidic solution directly to the pre-filled reagent.

The pre-filled reagent contains guanidinium salts, which, when combined with bleach can form highly reactive compounds. If any of these buffers are spilled, clean immediately with a suitable laboratory detergent and water. If the spilled liquid contains potentially infectious agents, clean the affected area first with laboratory detergent and water, and then with sodium hypochlorite at a concentration of 1% (v/v). The **Viral DNA and RNA Extraction Kit** comes with the following warnings and precautions.

Name of Component		Hazard pictograms (CLP)	Classification under CLP:	H- and P-statements
REAG 1	Lysis Buffer Washing Buffer A Washing Buffer B		Acute toxicity (oral), Category 4 Skin corrosion/irritation, Category 2 Serious eye damage/eye irritation, Category 2	<p>Hazard statements (CLP)</p> <p>H302: Harmful if swallowed. H315: Causes skin irritation. H319: Causes serious eye irritation.</p> <p>Precautionary statements (CLP)</p> <p>P264 - Wash hands, forearms and face thoroughly after handling. P280 - Wear protective gloves/protective clothing/eye protection/face protection/hearing protection. P321 - Specific treatment (see supplemental first aid instruction on this label). P337+P313 - If eye irritation persists: Get medical advice/attention. P501 - Dispose of contents/container to hazardous or special waste collection point, in accordance with local, regional, national and/or international regulation.</p>
	Elution Buffer	None	None	None

Please see MSDS for more details.

Precautions for Safe Handling

Do not dispose of the preparations or the packaging waste in drains leading to the sewage system or in the drainage system for waste not produced by industrial processing/analysis waste.

Any material in contact with reagents should be treated as a biological contaminant and treated in accordance with relevant local regulations.

Reagent Storage and Handling

The **Viral DNA and RNA Extraction Kit** should be stored at room temperature in a cool, dry and well-ventilated area. All components of the kit can be adequately stored for up to 18 months.

The kit should be used in a well-ventilated area, keep away from the source of heat, sparks, open flames, and smoking.

To avoid evaporation, the pre-filled reagent should be used immediately upon opening and should not be placed open for a long period of time.

Avoid exposure to UV light (e.g., for decontamination), which may result in accelerated aging.

Sample Handling and Storage

Avoid foam inside or on the samples. Depending on the starting material, sample pre-treatment may be required. Samples should be stored at room temperature (15~25°C) before starting the experiment.

Samples should be used immediately after collection to extract nucleic acid or stored at 2~8°C for further experiment within 24 hours. For long-term storage, the samples should be placed at -20°C.

For detailed information on sample pretreatment, please refer to 2.1.3.

Operation Guide

1. Automated Extraction Process

Automatic nucleic acid extractors (Libex and GeneRotex 96) enable nucleic acid extraction by magnetic beads. They use magnetic rods to move the beads adsorbed with nucleic acid into different reagent wells and then rapidly and repeatedly stirs the liquid through a mixing sleeve to mix the liquid and magnetic beads thoroughly. After cell lysis, nucleic acid adsorption, washing, and elution, the high-purity nucleic acid is obtained. Automatic nucleic acid extractors are characterized by high automation, rapid extraction speed, stable results, and ease of operation.

The user needs to load samples and magnetic bead nucleic acid extraction reagents into the reaction consumables, the nucleic acid extractors are going to perform all nucleic acid extraction operations according to the experimental procedures. Please refer to the user manual provided with an instrument for operating instructions.

2. Operation Steps of Automated Extraction

2.1 Automatic Nucleic Acid Extractor (model: Libex)

2.1.1 Edit Experiment Program

The extraction procedure of Libex Nucleic Acid Extractor is as follows:

No.	Well	Name	Waiting (s)	Mixing (s)	Magnet (s)	Speed	Volume (μL)	Heating State	Temp (°C)
1	1	Lysis	0	240	45	7	800	Closed	0
2	2	Washing 1	0	60	30	7	700	Elution	90
3	3	Washing 2	0	60	30	7	800	Elution	90
4	6	Elution	60	300	30	7	80	Elution	90
5	2	Release bead	0	60	0	7	700	Closed	0

2.1.2 Reagent Preparation

96-deep well plate:

Open the kit, take out the REAG1,, slowly invert it several times to resuspend the magnetic beads, then remove the plastic package and gently shake the 96-well plate so that the reagent and magnetic beads are concentrated on the bottom of the 96-well plate (a 96-well plate horizontal centrifuge can also be used for centrifugation at 500 rpm for 1 min). carefully tear off the aluminum foil sealing film before use to avoid liquid splashing.

6 strip tube:

Open the kit, take out the REAG1, slowly invert it several times to resuspend the magnetic beads, then gently shake the 6 strip tube so that the reagent and magnetic beads are concentrated on the bottom of the tube. Put the reagent on the docking (Note the direction and make sure that the tube is placed at the lowest level), carefully tear off the aluminum foil sealing film before use to avoid splashing, which is shown in Figure 2.

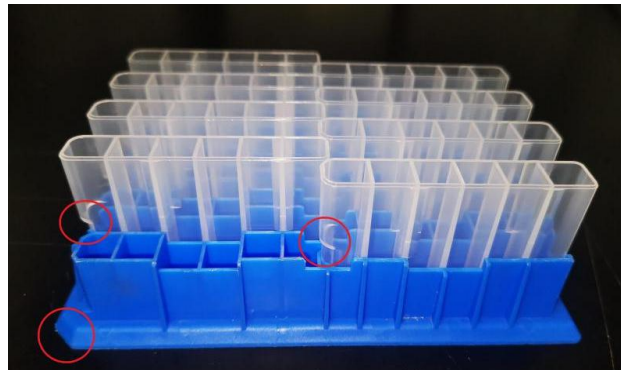


Figure 2. Put the 6 strip tube on the single kit docking

2.1.3 Adding Sample to the Reagent

96-deep well plate:

Add 200 μL of the sample that has been equilibrated to room temperature to column 1 or column 7 of the pre-filled reagent (Be aware of that column No. is for effective wells.)

6 strip tube:

Add 200 μL of the sample that has been equilibrated to room temperature to the first well of the pre-filled 6 strip tube.

⚠ Caution: When pipetting the sample, avoid having substance than liquid adhere to the tip of the sample injector; do not add the sample too quickly to avoid contaminating the upper portion of the well wall; and do not splash air bubbles to avoid contaminating adjacent wells.

⚠ Note: The following points should be taken into consideration when determining whether a sample is suitable for the Viral DNA and RNA Extraction Kit.

a. Type of sample: Swab samples etc.

b. Short-term storage: Samples can be used immediately after collection for nucleic acid extraction or stored at 2~8°C for testing with a maximum storage period of 24 hours.

c. Long-term storage: If the user does not operate the sample temporarily, it should be kept sealed in a refrigerator at -20°C.

2.1.4 Loading in Deep Well Plate

Place the 96-deep well plate or 6 strip tube in the Automatic Nucleic Acid Extractor and ensure the marked notch of the plate faces front.

Insert the mixing sleeve into the sleeve holder and close the cabin door.

⚠ Note: As shown in Figure 3 and Figure 4, the user should ensure that the 96-deep well plate and the single kit docking is properly positioned with the notch facing outward.

⚠ Note: Place the 96-deep well plate into the experiment cabin and push the magnetic rod covers into the right position. Check the position of the magnetic rod covers. Otherwise, instrument dysfunction or malfunction may occur and affect the experiment results.

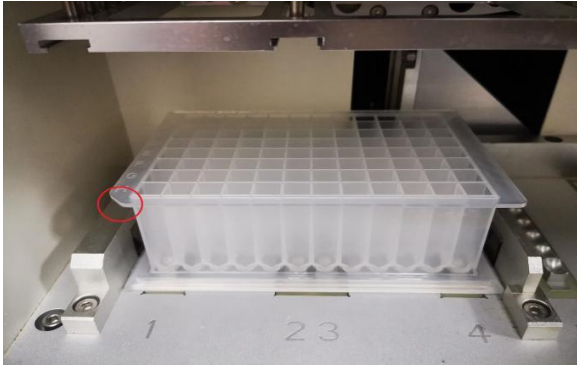


Figure 3. 96-deep well plate

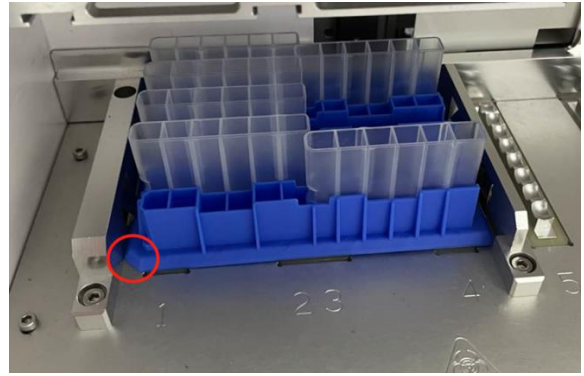


Figure 4. Put the single kit docking into the instrument

2.1.5 Experimental Procedure Run

For special operations please see 2.1.1. After the procedure is completed, the instrument will notify the user that the experiment has been completed. Transfer the extracted product from well 6 and well 12 to a clean centrifuge tube that is free of nuclease.

⚠ Note: If the user does not analyse the extracted product for the immediate use, please seal and store it in a refrigerator at -20°C.

⚠ Caution: Any used deep well plate and mixing sleeve should be considered as biological contaminants and disposed of in accordance with relevant regulations.

⚠ Caution: Using expired reagents or those that are not compatible with this instrument does not guarantee the expected results.

2.1.6 Cleaning and Maintenance of the Instrument

Follow the Cleaning and Maintenance of Instrument section in accordance with the instruction in the user manual provided with the equipment. Ensure that the experimental chamber is cleaned regularly to minimize the risk of cross-contamination.

2.2 Automatic Nucleic Acid Extractor (model: GeneRotex 96)

2.2.1 Edit Experiment Program

The extraction procedure of GeneRotex 96 Nucleic Acid Extractor is as follows:

Step	Name	Well	Stir (min:s)	Magnetic (min:s)	Wait (min:s)	Speed (rpm)	Volume (μL)	T Control (°C)
1	Lysis	1	04:00	00:45	00:00	2500	800	0
2	Washing 1	2	01:00	00:30	00:00	3000	700	100
3	Washing 2	3	01:00	00:30	01:00	3000	800	100
4	Elution	6	05:00	00:30	00:00	2500	80	100

2.2.2 Reagent Preparation

96-deep well plate:

Open the kit and take out the REAG1, slowly invert it several times to resuspend the magnetic beads. Gently shake the 96-well plate so that the reagent and magnetic beads are concentrated on the bottom of the 96-well plate (A 96-well plate horizontal centrifuge can also be used for centrifugation at 500 rpm for 1 min). Carefully tear off the aluminum foil sealing film to avoid the liquid splashing.

6 strip tube:

Open the kit and take out the REAG1, slowly invert it several times to resuspend the magnetic beads. Gently shake the 6 strip tube so that the reagent and magnetic beads are concentrated on the bottom of the tube. Put the reagent on the docking (Note the direction and make sure that the tube is placed at the lowest level), carefully tear off the aluminum foil sealing film before use to avoid the liquid splashing, which is shown in Figure 5.

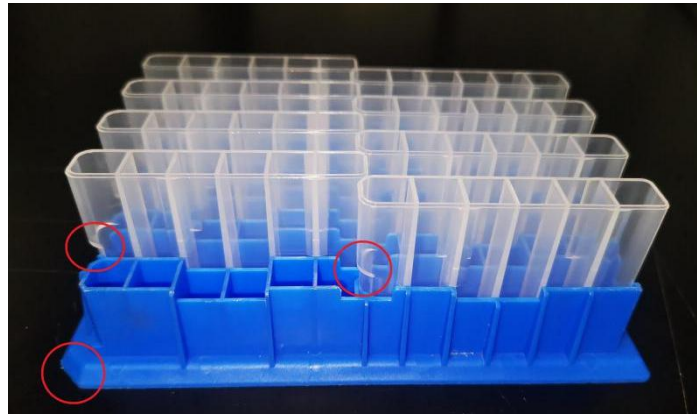


Figure 5. Put the 6 strip tube on the single kit docking

2.2.3 Adding Sample to the Reagent

96-deep well plate: Add 200 μ L of the sample that has been equilibrated to room temperature to column 1 or column 7 of the pre-filled reagent (Note the column No. is for effective wells.)

6 strip tube: Add 200 μ L of the sample that has been equilibrated to room temperature to the first well of the pre-filled 6 strip tube.

⚠ Caution: When pipetting the sample, avoid having substance than liquid adhere to the tip of the sample injector; do not add the sample too quickly to avoid contaminating the upper portion of the well wall; and do not splash air bubbles to avoid contaminating adjacent wells.

⚠ Note: The following points should be taken into consideration when determining whether a sample is suitable for the Viral DNA and RNA Extraction Kit.

a. Type of sample: Swab samples etc.

b. Short-term storage: Samples can be used immediately after collection for nucleic acid extraction or stored at 2~8°C for testing with a maximum storage period of 24 hours.

c. Long-term storage: If the user does not operate the sample temporarily, it should be kept sealed in a refrigerator at -20°C.

2.2.4 Loading in the Deep Well Plate

Properly position the 96-deep well plate or 6 strip tube containing the sample in the experimental cabin of the fully automatic nucleic acid extractor (GeneRotex 96).

⚠ Note: The user should ensure the marked notch of the plate of the 96-deep well plate and 6 strip tube are on the left, which is shown in Figure 6 and Figure 7.

Insert the rotatory mixing sleeve into column 2 and/or column 8 of the deep well plate and close the experimental cabin door.

⚠ Caution: The user must ensure that the rotary mixing sleeves are placed properly; otherwise, the instrument may operate abnormally, or the magnetic rods may become contaminated.

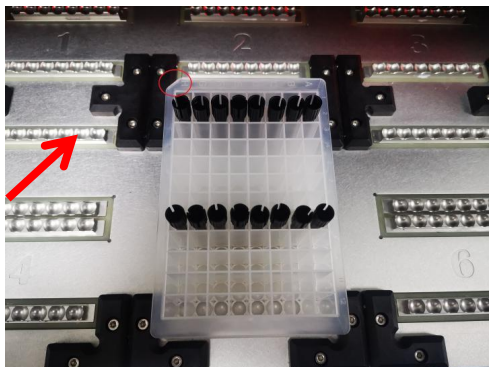


Figure 6. 96-deep well plate

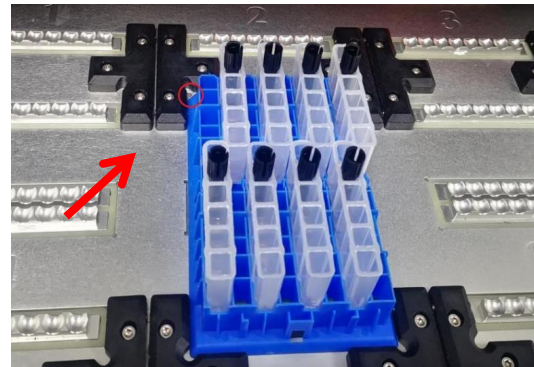


Figure 7. Put the single kit docking into the instrument

2.2.5 Experimental Procedure Run

For special operations please see 2.2.1. After the procedure is completed, the instrument will notify the user that the experiment has been completed. Transfer the extracted product from well 6 and well 12 to a clean centrifuge tube that is free of nuclease.

⚠ Note: If the user is not going to use the extracted product immediately, please seal and store it in a refrigerator at -20°C.

2.2.6 Cleaning and Maintenance of the Instrument

Follow the Cleaning and Maintenance of Instrument in accordance with the user manual provided with the equipment. Ensure that the experimental cabin is cleaned regularly to minimize the risk of cross-contamination.

Troubleshooting Guide

This troubleshooting guide should assist you in resolving any problems that arise during the experimental process. For more information and Frequently Asked Questions, please visit our Technical Support Center at <http://www.medtl.net>. The scientists in our Tianlong company's Technical Services Department are always available to answer any questions you may have about the information and protocols contained in the manual, sample and assay technologies (Contact information is included on the back cover or at <http://www.medtl.net>).

When an exception or error occurs during the experiment, the current run step is terminated/stopped. After resolving the error or exception, restart the run from the beginning. The troubleshooting guide is shown in the following table.

No.	Fault Symptom	Fault Cause	Handling Method
1	The well plate vibrates and the liquid splashes when tearing off the aluminum foil sealing film.	When tearing the film, please press the well plate to prevent it from rocking.	The reagent for this plate shall be scrapped, and re-extraction shall be performed.
2	Add the sample to unexpected wells.	Please read this manual carefully before adding samples.	The reagent for this plate shall be scrapped, and re-extraction shall be performed.
3	The amount of liquid in the reagent wells is insufficient.	/	Contact the after-sales service of our company.
4	Reuse of pre-filled components	Please read the precautions in this manual before using the kit.	Perform re-extraction of nucleic acid.
5	Abnormal noise from the instrument during extraction	The 96-deep well plate may be placed incorrectly.	Conduct reposition of the deep well plate.
		The rotary mixing sleeves may be inserted in wrong place.	Reinsert the mixing sleeve.
6	Poor extraction performance	Please follow the operation guide in the manual.	Contact the after-sales service of our company.
		The temperature control components of the instrument may be abnormal.	Contact the after-sales service of our company.

		Other	Contact the after-sales service of our company.
--	--	-------	---

* Ensure that the reagents have been preserved and used according to the manufacturer's instructions.

Quality Control

In accordance with Tianlong Company's ISO-certified Quality Management, each lot of **Viral DNA and RNA Extraction Kit** is tested against predetermined specifications to ensure consistent product quality.

Limitations of Test Methods














The system performance has been established through performance evaluation studies using swab samples to purify viral DNA and RNA.




It is users' responsibility to validate system performance for any procedures used in their laboratory that are not covered by the performance evaluation studies of Xi'an Tianlong Science and Technology Co., Ltd.

Although the kit is intended for use in public health and scientific research, the purity and quality of extraction results are also affected by the testing instruments and personnel. Moreover, the kit uses a specially formulated eluent that can affect the absorbance value, so it is not recommended to use a UV rays' spectrophotometer to measure the extraction effect directly.

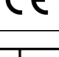The extraction kit is intended for use with clinical diagnostic samples, forensic materials, and scientific research samples. The instrument and operator have an effect on the concentration and purity of the extracted product. Any generated diagnostic results must be interpreted in conjunction with the other clinical or laboratory findings.

Safety Symbols and Signs

No.	Symbol	Implication
1		Catalogue number
2		Batch code
3		Contains sufficient for <N> tests
4		Use by date
5		Caution
6		Temperature limit
7		In vitro diagnostic medical device
8		Reminder
9		Manufacturer
10		Do not re-use
11		Conformed with EU standard
12		Authorized representative in the European Community
13		Content of the kit

14		Pre-filled 96-deep well plate/6 strip tube
15		Warning
16		PAP21: Not-corrugated cardboard

Contact Information

For technical assistance and more information, please contact our Technical Support Center at +86-29-82682132 (Tel), +86-29-82216680 (Fax), inquiry@medtl.com or contact your local distributor.

For a patient/user/third party in the European Union and in countries with similar regulatory regime (Regulation 2017/746/EU on IVD Medical Devices); if, during the use of this device or as a result of its use, a serious incident has occurred, please report it to the manufacturer and/or its authorised representative and to your national regulatory authority.

For up-to-date licensing information or product-specific disclaimers, please see the respective User Guide. Tianlong User Guides are available at www.medtl.net or can be requested from Tianlong Technical Services or the local distributor.

IFU_T014H/T016H/T060H_EN © 2023 Xi'an Tianlong Science and Technology Co., Ltd., all rights reserved.